Silk'n®

LED NECK MASK

FAC02


- EN USER MANUAL
- DE BEDIENUNGSANLEITUNG
- FR MODE D'EMPLOI
- NL GEBRUIKSAANWIJZING
- ES MANUAL DE USUARIO
- IT MANUALE D'USO
- TR KULLANIM KILAVUZU
- PT MANUAL DO UTILIZADOR

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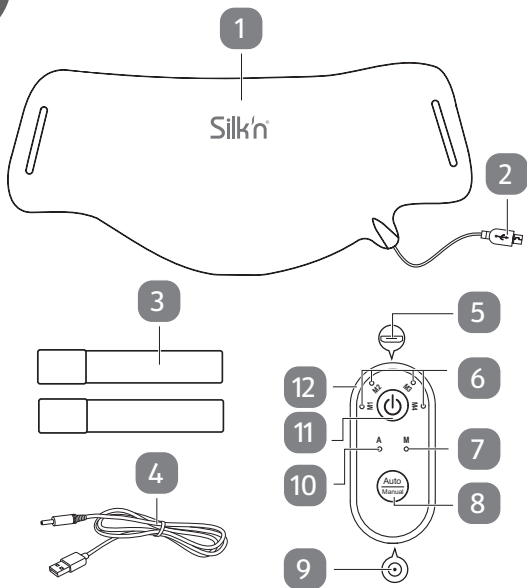
B

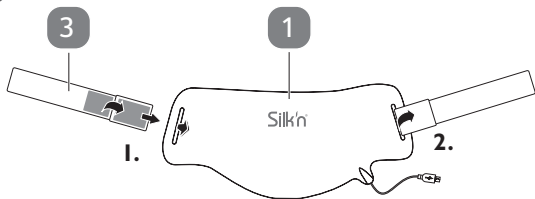Legal notice

Invention Works B.V. reserves the right to make changes to its products or specifications to improve performance, reliability, or manufacturability. Information furnished by Invention Works B.V. is believed to be accurate and reliable at the time of publication. However, Invention Works B.V. assumes no responsibility for its use. No license is granted by its implication or otherwise under any patent or patent rights of Invention Works B.V.

No part of this document may be reproduced or transmitted in any form or by any means, electronic or mechanical, for any purpose, without the express written permission of Invention Works B.V. Data is subject to change without notification.

Silk'n Beauty Ltd. has patents and pending patent applications, trademarks, copyrights, or other intellectual property rights covering subject matter in this document. The furnishing of this document does not give you any license to these patents, trademarks, copyrights, or other intellectual property rights except as expressly provided in any written agreement from Silk'n Beauty Ltd. Specifications are subject to change without notice.

Silk'n and the Silk'n logo are registered trademarks of Silk'n Beauty Ltd.

Invention Works B.V.
Galileïstraat 17, 3029 AL Rotterdam, The Netherlands
www.silkn.eu – info@silkn.eu

Table of contents

I	Package contents/device parts _____	5
2	General information _____	5
2.1	Reading and storing the user manual _____	5
2.2	Explanation of symbols _____	5
3	Safety _____	6
3.1	Proper use _____	6
3.2	Safety instructions _____	6
4	Checking the product and package contents _	10
5	Charging _____	10
6	Assembly _____	11
7	Operation _____	11
7.1	Carrying out the session _____	11
7.2	Manual mode _____	12
7.3	Auto mode _____	13
8	Regimen _____	13
9	Cleaning _____	14
10	Storage _____	14
11	Troubleshooting _____	14
12	Technical data _____	15
13	Disposal _____	15
13.1	Disposing of the packaging _____	15
13.2	Disposing of the product _____	15
14	Warranty _____	16
15	Customer service _____	16

I Package contents/device parts

- | | |
|-------------------------------|-------------------------------|
| 1 Neck mask | 7 Manual mode indicator light |
| 2 Plug | 8 Mode button |
| 3 Elastic strap, 2× | 9 Charging port |
| 4 USB cable | 10 Auto mode indicator light |
| 5 Connector | 11 Operation button (⏻) |
| 6 Program indicator light, 4× | 12 Remote control |

2 General information

2.1 Reading and storing the user manual

This user manual accompanies this LED Neck Mask (hereafter referred to as the "product"), and contains important information on setup and handling.

Before using the product, read the user manual carefully. This particularly applies to the safety instructions. Failure to do so may result in personal injury or damage to the product. The user manual is based on the standards and rules in force in the European Union and United Kingdom. When abroad, you must also observe country-specific guidelines and laws. Store the user manual for further use. Make sure to include this user manual when passing the product on to third parties.

2.2 Explanation of symbols

The following symbols and signal words are used in this user manual, on the product or packaging.

WARNING!

This signal symbol/word designates a hazard with moderate risk, which may result in death or severe injury if not avoided.

NOTICE!

This signal word warns of possible damage to property.

This signal symbol provides you with useful additional information on handling and use.

Read the user manual.

Declaration of Conformity: Products labelled with this symbol meet all applicable provisions of the European Economic Area.

Declaration of Conformity: Products labelled with this symbol meet all applicable provisions of the United Kingdom.

This symbol identifies voltage and the current rating for direct current.

Devices marked with this symbol are only suitable for household use (dry indoor rooms).

This symbol indicates the polarity of the charging port.

3 Safety

3.1 Proper use

The product is exclusively designed to recover your skin's glow, elasticity and texture. It also helps to reduce the appearance of fine lines, wrinkles and other skin imperfections. It is only intended for private use and not suitable for commercial purposes.

Only use the product as described in this user manual. Any other use is considered improper and may result in damage to property or persons. The manufacturer cannot be held liable for damages or injury incurred through improper or incorrect use.

3.2 Safety instructions

WARNING! RISK OF ELECTRIC SHOCK!

- Only use the product with the supplied USB cable.
- Only connect the USB cable of the product to a suitable mains adaptor (IEC/EN 60335)

certified SELV power supply) connected to an easily accessible socket so that you can quickly disconnect it from the mains supply in the event of a problem.

- Do not use the product if it is damaged or if the USB cable is defective.
- Do not open the housing; instead, have a qualified professional perform repairs. Contact the customer service team for this. Liability and warranty claims are waived in the event of repairs performed by the user, improper connection or incorrect operation.
- Only parts that comply with the original device data may be used for repairs. This product contains electrical and mechanical parts which are essential for providing protection against sources of danger.
- Do not immerse the product or USB cable in water or other liquids. Keep them away from water.
- Always switch the product off when you are not using it, when you intend to clean it, or in the event of a malfunction. Always unplug the USB cable when you are not charging the product.
- Never touch the USB cable with wet or damp hands.
- Never insert any objects into the housing.

- Keep all parts of the product away from open flames and hot surfaces.
- Lay the USB cable so that it does not pose a tripping hazard.
- Do not pull the USB cable out of the mains adaptor by the cord; instead, always pull it out by the cable plug.
- Do not kink the USB cable and do not lay it over sharp edges.
- Only use the product indoors.
- Do not use the product while charging.

 **WARNING! RISK OF INJURY!**

- The product can be used by persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning the use of the product in a safe way and understand the hazards involved. Children must not play with the product. Cleaning and user maintenance must not be performed by children.
- Keep children away from the product and USB cable.
- The product cannot be used by pregnant women, children, people with epilepsy, thyroid patients, people using thyroid medications, people who are photoallergic or use light allergy medication.

- The product cannot be used on skin diseases, wound ulceration, trauma and post-operative wounds.
- Do not use the product for more than 15 minutes per session.
- Consult your physician if you are uncertain whether you can use the product.
- Stop using the product immediately if you feel any discomfort, develop any abnormal reactions (such as itching or uncomfortable inflammation) or in case your skin becomes warm to the point of discomfort.
- Under abusive conditions, liquid may be ejected from the battery; avoid contact. Liquid ejected from the battery may cause irritation or burns.
- If battery fluid comes into contact with skin, eyes or mucous membrane, wash and rinse the affected areas immediately with clean water. Consult a physician if necessary.
- The battery must not be disassembled, short-circuited or thrown into a fire.

NOTICE! RISK OF DAMAGE!

- Keep the product away from open fire, hot surfaces, high temperatures or the effects of the weather (e.g. rain).
- Stop using the product if its parts exhibit cracks or breaks or are deformed. Only

replace damaged parts with corresponding original spare parts.

4 Checking the product and package contents

NOTICE! RISK OF DAMAGE!

- If you are not cautious when opening the packaging with a sharp knife or other pointed objects, you may quickly damage the product. Be very careful when opening it.
1. Take the product out of the packaging.
 2. Check to make sure that all parts are included (see **Fig. A**).
 3. Check whether the product or individual parts are damaged. If this is the case, do not use the product. Contact our customer service.

5 Charging

WARNING! RISK OF INJURY!

- Charge only with the supplied USB cable.
- Only charge the product at temperatures between 0 and 40 °C.
- Always unplug the product after charging has been completed.

- Fully charge the product before first use and when the battery power is low.
- When the battery power is low, the indicator lights on the remote control flash 10 times before switching off automatically. Recharge the product.
- The product cannot be used during charging.
- If the product is to be stored for an extended period, the battery level should be checked and the product should be charged regularly.

1. Connect the plug **2** of the neck mask **1** with the connector **5** of the remote control **12**.
2. Plug the smaller end of the USB cable **4** into the charging port **9** of the remote control.
3. Plug the larger end of the USB cable into the USB port of a compatible USB power source (e.g. USB power adaptor).
4. Plug the USB power adaptor into a mains socket for charging. Charging starts automatically. Check the program indicator lights **6** for the charging progress:
One indicator light = 25 % charged
Two indicator lights = 50 % charged
Three indicator lights = 75 % charged
Four indicator lights = fully charged
5. When charging is done, unplug the USB power adaptor from the mains socket and disconnect the USB cable from the charging port of the remote control.

6 Assembly

1. Open the hook-and-loop fasteners on the elastic straps **3** and guide the elastic straps through the side openings of the neck mask **1** (see step 1 in **Fig. B**).
2. Close the hook-and-loop fasteners to secure the elastic straps on the neck mask (see step 2 in **Fig. B**).

7 Operation

WARNING! RISK OF INJURY!

- Do not use the product for more than 15 minutes per session.

7.1 Carrying out the session

1. Fully charge the product and properly attach the elastic straps **3** (see chapter “Charging” and “Assembly”).
2. Thoroughly clean and dry your neck. Make sure it is free from any make-up or dirt.
3. Sit or lie down and place the neck mask **1** over your neck.
4. Use the hook-and-loop fasteners on the ends of the elastic straps to fit the neck mask comfortably on your neck.

- Connect the plug **2** of the neck mask into the connector **5** of the remote control **12**.
- Press and hold the operation button (⏻) **11** to switch on the product.
By default, the product goes into Manual mode with the M1 (red) program selected after switching on.
- If desired, press the mode button **8** repeatedly to select between the Auto mode and Manual mode.
The manual mode indicator light **7** or auto mode indicator light **10** turns on according to your selection.
To use the Manual mode, see chapter “Manual mode”.
To use the Auto mode, see chapter “Auto mode”.
- After the session, remove the neck mask and apply a serum and/or day/night cream on your neck.

7.2 Manual mode

- By default, the product goes into Manual mode with the M1 (red) program selected after switching on.
- If another program (M2, M3 or M4) is desired, press the operation button (⏻) **11** 1–3 times for selection. The corresponding program indicator light **6** turns on to show your selection.
Note: The product switches off when the operation button (⏻) is pressed a fourth time.
- When a program is selected, the manual mode indicator light **7** flashes 5 times before the indicator light stays lit and the product starts the selected program.
- Functions of the programs:

Program	LED	Best for
M1	Red	Reduces fine lines and wrinkles, stimulates collagen and helps to improve blood circulation
M2	Blue	Reduces acne-related redness and imperfections for clearer skin
M3	Yellow	Improves skin hydration and soothes dry and sensitive skin

Program	LED	Best for
M4	Purple	Anti-ageing and reduction of skin imperfections to improve the skin's vitality and radiance

- The session lasts for 10 minutes. The product switches off automatically after the session time has elapsed.
- You can also press and hold the operation button (⏻) to switch off the product any time.

7.3 Auto mode

- When Auto mode is selected, the auto mode indicator light **10** flashes 5 times before the indicator light stays lit and the product starts the session.
- In Auto mode, the product carries out the programs in below order:
 - M1 (red): 5 minutes
 - M2 (blue): 3 minutes
 - M3 (yellow): 4 minutes
 - M4 (purple): 3 minutes
- The whole session lasts for 15 minutes. The product switches off automatically after the session time has elapsed.
- You can also press and hold the operation button (⏻) **11** to switch off the product any time.

8 Regimen

If you are using the product for the first time, follow the schedule below:

First week	Use once every 3 days, max. 10 minutes at a time
Second week	Use once every 2 days, max. 10 minutes at a time
Third week	Use once every 2 days, max. 10–15 minutes at a time
Fourth week	Use once every day, max. 10–15 minutes at a time

As soon as results show, adjust the frequency to once or twice per week.

9 Cleaning

WARNING! RISK OF ELECTRIC SHOCK!

- Always unplug the product from the power source before cleaning.

NOTICE! RISK OF DAMAGE!

- Do not use any aggressive cleaners, brushes with metal bristles, or sharp or metallic cleaning utensils such as knives, hard scrapers, and the like.
- Do not use hot water for cleaning.
- Use a damp cloth or a water-based wipe to wipe off the neck mask . Then wipe it dry with a dry cloth.

10 Storage

All parts must be completely dry before being stored.

- Store the product in a clean and dry place, avoiding heat and direct sunlight.
- Store the product so that it is not accessible to children.
- For best protection, store the product in its original packaging.

11 Troubleshooting

Some problems may be caused by minor faults that you can fix yourself. To do so, follow the instructions in the following table. If it is still not possible to resolve the problem, contact customer service. Do not repair the product yourself.

Fault	Possible cause and solution
The product does not operate.	<ul style="list-style-type: none"> • Make sure the product is properly charged. • Make sure the product is not connected to a USB power source. The product does not work while charging.

Fault	Possible cause and solution
The product does not charge.	<ul style="list-style-type: none"> Make sure the USB cable 4 is securely connected between the charging port 9 and the USB power source.

12 Technical data

Model:	FAC02
Rated input:	5 V === 1 A
Working voltage:	3.7 V
Battery capacity:	1,500 mAh, 5.55 Wh (lithium)
Number of LED lights:	157
Weight:	169 g
Dimensions (L × W × H):	237.5 × 108 × 7.1 mm

13 Disposal

13.1 Disposing of the packaging

Sort the packaging before you dispose of it. Dispose of paperboard and cardboard with the recycled paper service and wrappings with the appropriate collection service.

13.2 Disposing of the product

(Applicable in the European Union and other European countries with separate collection systems of recyclable materials)

Old appliances may not be disposed of in the household waste!

Should the product no longer be capable of being used at some point in time, **dispose of it in accordance with the regulations in force in your city or country.** This ensures that old appliances are recycled

in a professional manner and also rules out negative consequences for the environment. For this reason, electrical equipment is marked with the symbol shown here.

Batteries and rechargeable batteries may not be disposed of with household waste!

As the end user you are required by law to bring all batteries and rechargeable batteries, regardless whether they contain harmful substances* or not, to a collection point run by the communal authority or to a retailer, so that they can be disposed of in an environmentally friendly manner.

Bring the entire product (including the rechargeable battery) to your collection point and make sure that it is in an uncharged state!

* labelled with: Cd = cadmium, Hg = mercury, Pb = lead

14 Warranty

This product is covered by a 2-year warranty according to European regulations and law. The extent of the warranty on this product is limited to technical defects caused by faulty production processes. In the event you want to claim warranty please be sure to contact our customer service for instructions. They may be able to solve your problem without needing to return the product to the store or our service center. Our customer service is always pleased to help you!

15 Customer service

For more information about the products of Silk'n visit your regional Silk'n website: www.silkn.eu. If the device is damaged, defective, in need of repair or you need our assistance, please contact your nearest Silk'n service center. This manual is also available as a PDF download from: www.silkn.eu.

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© 2023 Silk'n Beauty Ltd. All rights reserved. Silk'n and the Silk'n logo are registered trademarks of Silk'n Beauty Ltd.
Silk'n UK Ltd.