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ISH702

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**Dittmann.**  
HEALTH

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SHENZHEN DONGDIXIN TECHNOLOGY CO.,LTD., FLOOR 1-2,  
NO.3 BUILDING, FANSHEN XUSHENG INDUSTRIAL ESTATE  
XILIXIAOBAIMANG 518108 NANSHAN DISTRICT,  
SHENZHEN P. R. CHINA

Bitte Gonen  
Modell: CH13000  
Type:Pro

MEDNET SWISS GMBH  
D4 PLATZ 4, 6039 ROOT D4, SWITZERLAND  
MEDEVVOY SWITZERLAND,  
GOTTHARDSTRASSE 28, 6302 ZUG, SWITZERLAND

SHANGHAI INTERNATIONAL HOLDING CORP. GMBH (EUROPE)  
EIFFESTRASSE 80, 20537 HAMBURG, GERMANY

HANDELSHAUS DITTMANN GMBH  
KISSINGER STRASSE 68  
97727 FUCHSSTADT, GERMANY

ISH 702

**Dittmann.**

**CE** 0598

Copyright © Dittmann GmbH, 2024

Dok./REV.-Nr. SM431\_20240114

# **INSCET BITE HEALER**

**ISH702**

*Relieve itching and  
swelling quickly*

**Dittmann.**  
HEALERS

EN

# TABLE OF CONTENTS

## TABLE OF CONTENTS

Introduction .....	3
Important safety precautions and warnings .....	5
Know your device .....	11
Inserting batteries.....	13
Easy steps to get started .....	14
Other important functions.....	17
Specifications .....	18
Cleaning and storage .....	19
Disposal .....	20
Troubleshooting .....	21
Glossary of symbols .....	22
Important information regarding electromagnetic compatibility (EMC) .....	24
Warranty .....	25

# INTRODUCTION

**EN**

## INTRODUCTION

Thank you for purchasing the for relieving itching and swelling of insect bites.

In order to use the device safely, read the complete manual carefully before using the device for the first time.

Keep this instruction manual in a convenient place, or store with the device for future reference.

The device is a medical device designed to be used in the home healthcare environment.

It is suitable for use by all who can control the device and understand the instructions.

Patient groups: 3 years old and older

Intended users: Lay users, for children under the age 12 should be use under adult supervision.

### **operating principle**

Bite Gone can relieve the insect bite/sting skin by concentrating heat on the affected area. The mechanism of action of applying concentrated heat (local hyperthermia) is based upon the brief,

**EN**

## INTRODUCTION

concentrated effect of heat on a very small area of skin. Even a locally limited application of heat may be enough to trigger a pain and itch-relieving reaction of the body. This means that once the device heats to the right temperature and you use it on the area of skin affected by the sting or bite, that heat transfer quickly reduces skin itching, pain and swelling.

### Indications for use

The device is intended for symptomatic treatment of itching, pain and swelling caused by insect stings and bites.

### Standard parts

The package contains the following components:

- 1x Insect bite healer
- 1x instructions for use
- 2x batteries, type AAA (LR03 1.5V)

# IMPORTANT SAFETY PRECAUTIONS AND WARNINGS

**EN**

## IMPORTANT SAFETY PRECAUTIONS AND WARNINGS

It is important that you read all the contraindications, warning and precautions included in this manual because they are intended to keep you safe, prevent injury and avoid a situation that could result in damage to the device.

### SAFETY SYMBOLS USED IN THIS MANUAL

#### **WARNING**

Indicates a potentially hazardous situation which, if not avoided, could result in serious injury and equipment damage.

#### **CAUTION**

Indicates a potentially hazardous situation which, if not avoided, may result in minor or moderate injury to the user, or damage to the device or other property.

For a patient/user/third party in the European Union and in countries with identical regulatory regime (Regulation (EU) 2017/745 on Medical Devices); if, during the use of this device or as a result of its use, a serious incident has occurred, please report it to the manufacturer and/or its authorized representative and to your national authority.

**EN**

## IMPORTANT SAFETY PRECAUTIONS AND WARNINGS

### Contraindications

- Do not use if you suffer from diabetes, The pain may have diminished, and diabetics may burn themselves.
- Not for use on children under the age of 3.
- Do not use on inflamed skin, reddened skin, irritated skin, rashes, burns, frostbite, bruises, swellings, open or healing wounds, unhealed or fresh post-operative scars.
- Do not use, if suffering with fever.
- Do not use the device on sensitive areas of skin.
- Do not use, if you have lotions, creams, and gels on your skin at the time of insect bite.
- Do not use, if your perception of pain is reduced (i.e. metabolic disorders).
- Do not use if you suffer from persistently inflamed skin due to long term health application to the same skin area (i.e. where the insect bite occurred).

### General warning

- The patient may be the proper operator of the device.
- Not for use on children under the age 12 without adult supervision.
- Not for use by disabled people and people with reduced sensitivity to heat, such as people with illness related skin alterations and scars in application area.

## IMPORTANT SAFETY PRECAUTIONS AND WARNINGS

**EN**

- Do not use if you have thromboembolism or recurrent malignant growth or unexplained chronic pain in any area of the body.
- Do not use near the eyes, eyelids, mouth, or mucus membranes.
- Ensure that you observe general rules of hygiene, when using the Bite Gone.
- Seek medical advice if bite symptoms persist or worsen despite using the Bite Gone.
- Seek medical advice immediately, if the bite is a tick bite as it may result in a transmission of pathogens (i.e. tick-borne encephalitis (TBE) or borreliosis).
- Seek medical advice is the insect bite related symptoms, such itching increases and skin redness, swelling, dizziness, nausea, shortness of breath, formation of wheals, or drop or/increase in blood pressure occurs.
- If in doubt as to your suitability for using the Bite Gone, consult a medical practitioner (i.e. a nurse or a doctor).
- Make sure children do not treat the device as a toy.
- Keep equipment out of reach of children and pets.
- Children principle Packaging Materials (choking risk).
- Observe general hygiene practices when using the equipment.
- This device is not a substitute for medical consultation and treatment. So consult your doctor if you have any pain or illness.
- If the device does not work properly, if you feel uncomfortable or painful, please stop using it immediately.

**EN**

## IMPORTANT SAFETY PRECAUTIONS AND WARNINGS

- Do not service and maintain while the device is in use.
- Do not modify the device without authorization of the manufacturer. This could cause improper functioning. You can replace battery and clean the device, other maintenance is not allowed.
- If you have allergic reaction to materials, please stop treatment immediately and consult your physician.
- Please keep the battery compartment closed during operation and avoid touching the battery and the body part at the same time.

### Matters needing attention

- Do not touch the device with wet hands. Do not inject water into the device. The equipment can only be operated in a completely dry state.
- Do not use detergent or solvent to clean or maintain the device. The injected liquid can damage electronics and batteries and may cause device failure.
- The device can only be used by one person. Multiple users are not recommended.
- Protect the device from direct sunlight and store it dry.
- Protect equipment from shock, moisture, dust, chemicals, severe humidity fluctuations, and near heat sources (furnaces, radiators).

## IMPORTANT SAFETY PRECAUTIONS AND WARNINGS

**EN**

- Make sure no water enters the equipment. If this happens, remove the battery and use it only when the device is completely dry.
- Do not place heavy objects on the storage device.
- Make sure there is no visible damage before use. If in doubt, do not use, please contact your dealer or customer service address.
- Equipment must be maintained in expected environmental conditions.
- Do not leave the device unattended during use. If something goes wrong, stop using it immediately.

### Battery precautions

#### Warning

- If battery leakage comes into contact with skin or eyes, wash the affected area with water and seek medical advice.
- Please keep the battery away from children. There is a risk of suffocation if swallowed.
- Keep away from the source of fire, which could cause an explosion.
- When the battery runs out and leaks, wear protective gloves and clean the battery compartment with a dry cloth.
- Do not remove, open or crush the battery.

**EN**

## IMPORTANT SAFETY PRECAUTIONS AND WARNINGS

### Caution

- Note the polarity plus (+) and minus (-) signs.
- Protect the battery from heat.
- Batteries are not allowed to be charged or short-circuited.
- If the device is not used for a long time, take the battery out of the battery compartment.
- Use only the same battery.
- Be sure to replace all batteries at the same time.

### Side effects

Depending on sensitivity, the treatment temperature may cause some discomfort. In sensitive individuals, a temporary reddening of the skin may occur after use. In rare cases, users with sensitive skin types may experience skin irritation.

### Residual risk

The heated ceramic surface may reach a maximum temperature of 127°F (53°C). You may get skin burn if used on the same skin area more than once within a short period of time.

## KNOW YOUR DEVICE

**EN**

### KNOW YOUR DEVICE

#### Features

- The mosquito bites and insect stings hardly avoided if you keep yourself outdoors for a long time in the summer.
- Bite Gone helps to relieve itching from mosquito bites and insect stings quickly and precisely and accelerates healing process.
- Thanks to its handy and compact design, the device can be taken anywhere - whether for hiking, camping or travelling.
- The treatment of mosquito bites and insect stings with the device is only done by heat and without any chemical substances.
- The application is therefore also suitable for whole family including pregnant women.
- One key operation.
- 3 different programs.
- Drug and chemical free.
- Relieve and accelerate mosquito and insects bite itching and swelling quickly.
- Fast heating, ceramic heating plate.
- Accelerate healing process.
- Handy design for on-the-go.
- User friendly, simple and easy operation.
- There are lanyard holder for easier carrying.


EN

## KNOW YOUR DEVICE

### Front and rear panel

- 1) Indicator
- 2) Function Button (Switch On/Off Button)
- 3) Torch light
- 4) The heating head
- 5) Battery compartment
- 6) Lanyard holder

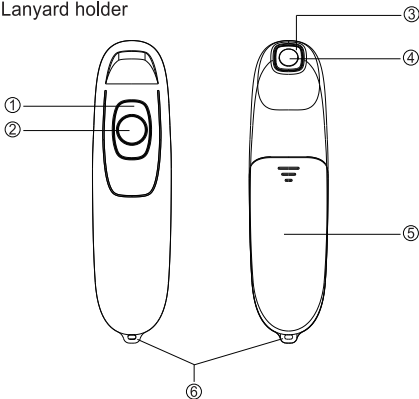
## INSERTING BATTERIES

EN

### INSERTING BATTERIES

1. Remove the battery cover on the back of the device as shown by the graphic to the right.
2. Insert 2xAAA batteries. Make sure the positive (+) and negative (-) signs correspond with the markings in the device when inserting batteries. Reinstall the battery cover as shown by the graphic to the right.
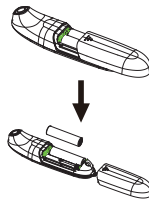
#### Notes:

- Please use 2xAAA batteries in this unit.
- Remove the batteries if the device is not in use for long periods of time.
- Do not mix old and new batteries or different types of batteries.
- Warning: If batteries leak and come into contact with the skin or eyes, wash immediately with large amounts of water.
- Batteries must be handled by an adult. Keep batteries out of the reach of children.
- Remove exhausted batteries from the unit.
- Dispose the used batteries safely according to local regulations.

 EN

## EASY STEPS TO GET STARTED

### EASY STEPS TO GET STARTED

1. Prior to applying the device, ensure that the sting area is clean and that all remains of the sting have been carefully removed. Ensure that the device is cleaned before and after each use.  
For best outcomes, use the device immediately after the sting or bite occur. However, even at the later stage after the sting, the device can still be helpful for symptomatic relief and feeling of greater comfort.
2. Ensure that the batteries are fitted in the battery compartment and the compartment lid is closed.
3. Long press the button(2) to power up – you will see the Indicator(1) light up green, and the Torch light (3) will also light up.

**Caution:**

The time required for the device to warm from the minimum storage temperature (-10°C) between use until the device is ready for use at ambient temperature 20°C: about 2 hours.

The time required for the device to cool from the maximum storage temperature (55°C) between use until the device is ready for use at ambient temperature (20 °C) : about 2 hours

## EASY STEPS TO GET STARTED

**EN**

4. Select an appropriate programme:
  - For sensitive skin or first use - select program 1(3 seconds /48°C) by pressing the button once (1 x) .The indicator(1) will light blue.
  - For normal skin or mild symptoms - select program 2 (3 seconds/51°C) by pressing the button twice (2 x).The indicator(1) will light white.
  - For normal skin and more severe symptoms-select program 3 (6 seconds/51°C) by pressing the button thrice(3 x) .The indicator(1) will light purple.

**Caution:**

When using the device for the first time, you may opt to use it on the Program 1 to enable you to test your tolerance of the heat treatment.

5. At the end of the application, the device sounds an acoustic signal, Indicator (1) lights up green, and Torch light (3) also lights up. Remove the hot plate from the skin as soon as the acoustic signal sounds. Heating of the hot plate automatically stops.
6. If you do not intend to use the device, Long press the button to the power off.

**Caution:**

Alternatively, you can re-use the device on another sting. However, please note that the maximum number of 5 applications per hour on the same treatment area must not be exceeded.

EN

## EASY STEPS TO GET STARTED

7. Ensure you clean the device after use.

**Caution:**

Ensure that you follow cleaning instructions outlined in user manual. The Inscet bite healer is a multiple use device. In the interest of hygiene and prevent potential cross contamination, the device must be cleaned before and after each use. Consult user manual for details and helpful advice. Ensure that the device is switched off, cool and clean prior to placing in storage. If you intend to store it unused for a long time, it is best to remove batteries from the compartment for the duration of storage.

## OTHER IMPORTANT FUNCTIONS

EN

### OTHER IMPORTANT FUNCTIONS

#### **Automatic shutdown**

If no operation is performed for more than one minute, the device automatically shuts down.

#### **Low electricity prompt**

If the voltage of the point battery is too low, three prompt tones will sound and the Indicator (1) will light up green and blink at a frequency of 1Hz, automatic shutdown after 3 seconds.

#### **Abnormal temperature protection**

The device cannot be used when the temperature of the device is lower than 5 °C or higher than 53 °C (or the temperature sensor is damaged).

When these abnormal conditions occur, the treatment light will blink to alert.

**EN**

## SPECIFICATIONS

### SPECIFICATIONS

The insect bite healer technical specifications:

- Product name: Bite Gone
- Model: CH3000
- Type: Pro
- Power supply: 2 x 1.5 V AAA LR03 batteries
- Maximum temperature: 51°C ±2°C
- Functions: 3 x settings (Program 1: 3 Seconds/48°C, Program 2: 3 Seconds/51°C, Program 3: 6 Seconds/51°C)
- Expected operating duration: circa 300 applications
- Operating conditions: 10°C to 40°C, RH 30% - 75%; 700hPa-1060hPa
- Storage conditions: -10°C to 55°C, RH 10% - 90% 700hPa-1060hPa
- Dimensions: 116,2 x 30 x 26,1mm +/-0,5mm
- Weight: about 50g (with batteries)
- IP Classification: IP22
- Applied part: Heating head
- Classification of ME Equipment: Internally powered/Continuous
- Service life of the device: 3 years
- Service life of the battery: New batteries will last for approx.300 times (3 seconds program)
- Software version: V1

# CLEANING AND STORAGE

**EN**

## CLEANING AND STORAGE

### Cleaning the device

- Please observe hygiene guidelines at all times.
- It is recommended to clean the device before and after each use.
- Please shut down the device and remove the batteries before cleaning.
- The device temperature must be restored to room temperature during cleaning.
- Use a soft towel (a small amount of alcohol or other medical cleanser is recommended) to wipe the heater and the casing.
- Do not let detergent into the machine, Irreversible damage may be caused to the device.
- Do not immerse the device in detergent, Irreversible damage may be caused to the device.
- Do not use chemicals (such as benzene, thinner).
- Do not clean in washing machine or other cleaning equipment.
- Do not use hard brushes to clean the equipment, which may cause scratches.
- If the device is used between different skin areas, clean the device with alcohol or other medical cleanser to prevent the spread of bacteria from one skin area to another.

**EN**

## DISPOSAL

### Storing the device

- Place the unit, and manual back into retail box. Store the box in a cool, dry place, 14°F~131°F (-10°C~55°C); 10%~90% relative humidity. Failing to do so may affect the performance of the device.
- Do not store in places that can be easily reached by children.
- Store the device in the dry, clean place. It is detrimental to leave the product exposed to any chemical solvent, lint, dust, direct sunshine or high temperature.
- When not in use for a long period, remove the batteries before storage to avoid liquid discharge from batteries.

### DISPOSAL

Used fully discharged batteries must be disposed of in a specially labeled collection container, at a toxic waste collection point or through an electrical retailer. You are under legal obligation to dispose of batteries correctly.

Please dispose of the device in accordance with local regulations.

Batteries and accumulators that are not firmly enclosed in the electrical or electronic appliance and can be removed without destroying them must be separated from the appliance before it is handed in at a collection point and disposed of in a designated disposal facility.

# TROUBLESHOOTING

**EN**

## TROUBLESHOOTING

If the unit does not operate after taking these measures, contact the nearest dealer.

PROBLEM	POSSIBLE CAUSES	POSSIBLE SOLUTION
The unit cannot power on	Are the batteries exhausted?	Replace the batteries.
Equipment not hot enough	Battery low voltage	Replace the batteries
The Indicator lights up green and blinks three times, and then the machine shuts down	Battery low voltage	Replace the batteries
Although the battery has been replaced, it still won't turn on	The wiring may be defective	Please contact customer service.

EN

## GLOSSARY OF SYMBOLS

### GLOSSARY OF SYMBOLS

Electrical devices are recyclable material and should not be disposed of with household waste after their useful life. Help us protect the environment and save resources and by taking this device to the appropriate collection point. Please contact the organization which is responsible for waste disposal in your area if you have any questions.

Caution

Type BF Applied Part

Please refer to instruction manual

**IP22**

The first number 2: Protected against solid foreign objects of 12,5 mm  $\Phi$  and greater. The second number: Protected against vertically falling water drops when enclosure tilted up to 15°. Vertically falling drops shall have no harmful effects when the enclosure is tilted at any angle up to 15° on either side of the vertical.

# GLOSSARY OF SYMBOLS

**EN**

Serial number.

The name and the address of the manufacturer

Transportation and storage temperature  
from -10 C to 55 C

Transportation and storage humidity limits  
from 10% to 90%

Transportation and storage atmospheric  
pressure limits from 700hPa to 1060hPa

Date of Manufacture:  
2024-01 (year/month)

The name and the address of the Authorized  
EC-representative in Europe

Medical Device

Unique device identifier

Batch code: V0624CH3000

CE mark

**EN****EMC**

## **IMPORTANT INFORMATION REGARDING ELECTROMAGNETIC COMPATIBILITY (EMC)**

Important information regarding electromagnetic compatibility  
This device should not be used adjacent to or stacked with other equipment and that if adjacent or stacked use is necessary, this device should be observed to verify normal operation in the configuration in which it will be used.

Use of accessories, detachable parts and materials other than those specified or provided by the manufacturer of this device could result in increased electromagnetic emissions or decreased electromagnetic immunity of this equipment and result in improper operation.

Portable RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 30 cm (12 inches) to any part of the device, including cables specified by the manufacturer. Otherwise, degradation of the performance of this equipment could result. When the operating environment is relatively dry, strong electromagnetic interference usually occurs.

At this time, the device may be affected as follows:

- the device stops output;
- the device turns off;
- the device restarts;

The above phenomenon does not affect the basic safety and essential performance of the device, and the user can use it according to the instruction. If you want to avoid the above phenomenon, please use it according to the environment specified in the manual.

The device is classified as a Group 1, Class B devices according to CISPR 11. The device is suitable for home healthcare environments and so on. The device meets the immunity levels for home healthcare environment. The device has been tested according to the EN 60601-1-2 standard by an accredited lab and found to be compliant for each emissions and immunity standard or test.

# WARRANTY

**EN**

## WARRANTY

The device has been developed and manufactured with great care.

The statutory warranty period is 24 months as of the purchase date for material and manufacturing defects of the product. Please keep the receipt as proof of the purchase of the device in order to be able to assert a possible warranty claim. The warranty does not include:

- damage due to incorrect usage
- defects that were already known to the customer when purchasing the product
- wear and tear parts
- damage due to unauthorised interventions and personal negligence on the customer's part

After the expiration of the warranty period you may send a defective device to the address specified below for repair. Repairs after the expiration of the warranty period are subject to charge. Do not hesitate to contact us if technical problems, questions and warranty claims concerning this device arise:

**NOTE:**

Please contact the service centre in the event of a complaint concerning the device!

If necessary, the service centre will arrange for the collection of the device.

Only PREPAID parcels are accepted by the service centre!  
 Handelshaus Dittmann GmbH  
 Abteilung Service-Center  
 Kissinger Straße 68  
 D-97727 Fuchsstadt/Germany  
 E-Mail: hotline@servicecenter.tv  
 Telephone hotline: +49 (0) 180-6000228 (€ 0.20 per call from a German  
 landline  
 number, max. € 0.60 per call from a German mobile phone networks)  
 www.dittmann-gmbh.com  
 Yours sincerely

SHENZHEN DONGDIXIN TECHNOLOGY CO.,LTD., FLOOR 1-2,  
 NO.3 BUILDING, FANSHEN XUSHENG INDUSTRIAL ESTATE  
 XILIXIAOBAIMANG 518108 NANSHAN DISTRICT,  
 SHENZHEN P. R. CHINA

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MEDENVOY SWITZERLAND,  
 GOTTHARDSTRASSE 28, 6302 ZUG, SWITZERLAND

SHANGHAI INTERNATIONAL HOLDING CORP. GMBH (EUROPE)  
 EIFFESTRASSE 80, 20537 HAMBURG, GERMANY

HANDELSHAUS DITTMANN GMBH  
 KISSINGER STRASSE 68  
 97727 FUCHSSTADT, GERMANY ISH 702

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