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LIST OF CONTENTS

Items Supplied	9
At A Glance	9
Symbols	10
Signal Words	10
Further Explanations	10
Intended Use	11
When May the Device Not Be Used?	11
Safety Notices	11
Use	13
Cleaning and Storage	14
Troubleshooting	14
Technical Data	15
Disposal	15

ITEMS SUPPLIED

- 1 x massage device
- 1 x mains adapter
- 1 x operating instructions

Check the scope of delivery for completeness and the components for transport damage. If you find any damage, do not use the product but contact our customer service department.

Remove any possible films, stickers or transport protection from the product.

Never remove the rating plate and any possible warnings!

AT A GLANCE

(Picture A)

- 1 2 x 4 massage balls, heating and red light element
- 2 Control Panel
- 3  On/Off button: When switched on the massage device starts the massage and changes direction approx. every minute. The red light and heat function is automatically switched on at the same time.
- 4  change of direction button: The massage balls massage either clockwise or anticlockwise. After approx. 1 minute the direction changes.
- 5  change of speed button: the speed can be set to three levels (fast, medium, slow)
- 6  button: Red light and heat can be switched off or on again
- 7 Hand loops
- 8 Connecting socket for the mains adapter

Dear Customer,

We are delighted that you have chosen our **VITALmaxx Shiatsu massage device**.

If you have any questions about the device and about spare parts/accessories, please contact the customer service department via our website: www.ds-group.de/kundenservice

Information About the Operating Instructions

Before using the device for the first time, please read through these operating instructions carefully and keep them for future reference and other users. They form an integral part of the device. The manufacturer and importer do not accept any liability if the information in these operating instructions is not complied with.

SYMBOLS


Danger symbols: These symbols indicate possible dangers. Read the associated safety notices carefully and follow them.

Supplementary information

Read the operating instructions before use!

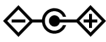Polarity of the coaxial connector (mains adapter)

Protection class II

Direct current

Alternating current

Energy efficiency class VI

Switching power supply

Safe electrical isolation

May only be used indoors!

SIGNAL WORDS

Signal words indicate hazards if the associated notices are not followed.

DANGER – high risk, results in serious injury or death

WARNING – medium risk, may have serious injury or death

CAUTION – low risk, may result in minor or moderate injury

NOTICE – may result in risk of damage to material

FURTHER EXPLANATIONS

Key numbers are shown as follows: (1)

Picture references are shown as follows: (Picture A)

INTENDED USE

- This device is intended to be used for loosening the neck and back muscles and delivers relaxation. The device may also be used to massage abdominal, buttock and leg muscles. It is not a replacement for medically required treatment.
- The device is for personal use only and is not intended for commercial applications.
- The device should only be used as described in these instructions. Any other use is deemed to be improper.
- The following are excluded from the warranty: all defects caused by improper handling, damage or unauthorised attempts at repair. This also applies to normal wear and tear.

WHEN MAY THE DEVICE NOT BE USED?

- The following people may use the device **only following consultation with a doctor**:
 - people who have no feeling in their neck, back, buttocks and / or thighs / lower legs
 - people who have back complaints or spinal disorders
 - people with circulation problems
 - women who are pregnant or have just given birth
 - people with a pacemaker
 - people who have undergone an operation for which the treatment is not yet complete.
- **Only use the device if you have no physical illnesses or ailments!** If you are unsure whether the device may be used, consult a doctor.
- Do not use the device near to or on injuries, swelling, skin irritations, infected skin, warts, fungal infections, abnormal changes in the skin or body parts with no feeling.

SAFETY NOTICES

WARNING: Note all safety notices, instructions, illustrations and technical data provided with this device.

Failure to follow the safety notices and instructions may result in electric shock, fire and / or injuries.

- This device should not be used by people (including children) with impaired sensory or mental faculties or who do not have adequate experience and / or knowledge of such devices, unless they are supervised by a person responsible for their safety, or have received instructions from such a person on how to use the device.
- Children must be supervised to ensure that they do not play with the device.
- Only ever use the device with the mains adapter that is supplied with it (LX-065-E120200). Do not use the mains adapter supplied with another device.

- When the heat function is switched on, the device has a hot surface. People who due to illness or for other reasons are not sensitive to heat must not use the heat function.
- The device may not be used by people who are vulnerable and are unable to react to overheating.
- Keep sharp and pointed objects away from the device to ensure that it is not damaged.
- The device must always be disconnected from the mains power when it is left unattended and before cleaning it.
- In the interest of your own safety, you should inspect the device for damage each time before you use it. It may only be used if it has no visible signs of damage!
- Do not make any modifications to the device. Do not replace the connection yourself. If the device or the connection are damaged, they must be replaced by the manufacturer, customer service department or a similarly qualified person (e. g. specialist workshop) in order to avoid any hazards.

DANGER – Danger of Electric Shock

- Use and store the device only in closed rooms.
- Do not operate or store the device in rooms with high humidity.
- Never immerse the device or the connection in water or other liquids and ensure that they cannot fall into water or become wet.
- If the device becomes wet or falls into water while it is connected to the mains power, pull the mains plug out of the plug socket immediately. **Only after this** should you retrieve the device from the water. Never attempt to pull it out of the water while it is connected to the mains power!
- Never touch the device or the connection with wet hands when the device is connected to the mains power.
- Do not use the device if it has malfunctioned or been dropped into water to avoid any hazards. Have it checked in a specialist workshop before using it again.

WARNING – Danger of Fire

- Do not use the device under a blanket or similar item or in the vicinity of combustible material.

WARNING – Danger of Injury

- Keep the packaging material away from children and animals. There is a danger of suffocation!
- Make sure that the connection is always kept out of the reach of small children and animals. There is a danger of strangulation.
- Do not insert any fingers or other body parts into the device. Do not open the zip during use.
- Make sure that the connected mains cable does not present a trip hazard. There is a danger of falling.

NOTICE – Risk of Damage to Material and Property

- In the event of damage / faults (e. g. unusual noises, odour or smoke development) and in a thunderstorm, switch the device off immediately and pull the mains adapter out of the plug socket.
- The device must be switched off when inserting or removing the mains adapter. Always pull the mains adapter and never the mains cable to disconnect the device from the mains power.
- Protect the device from heat, flames, direct sunlight, extreme temperatures, moisture and collisions.
- Keep the device in a dry place at temperatures of between 0 and 40 °C.

USE

Please Note!

- No pressure should be exerted on the device as this will slow down or stop the massage and the device could be damaged. Do not stand or lie on the device.
- Do not use the device for longer than 30 minutes without stopping.
- Do not massage one group of muscles for longer than 15 minutes.
- Switch off the device immediately if you experience pain or feel unwell during use.
- Only connect the device to a plug socket that is properly installed. The plug socket must also be readily accessible after connection so that the connection to the mains can quickly be isolated. The mains voltage must match the rating plate of the mains adapter.
- Only use suitable extension cables whose technical data is the same as that of the device.

- Make sure that the connecting cable is not squashed, bent or laid over sharp edges and does not come into contact with hot surfaces.
- Unwind the mains cable fully before connecting the mains adapter.
- For reasons of hygiene, avoid direct contact of the skin with the device wherever possible. It should not come into contact with skin that has had cream or oil applied to it.
- A tingling sensation/itching during or after the massage is normal due to the circulation-boosting effect.
- Do **not** fall asleep while you are using the device!

Do not use the device before going to bed. The massage has a stimulating effect, which may lead to problems getting to sleep. We recommend that you wait an hour after mealtimes before using the device.

Areas of Application

When this device performs a massage, four massage heads massage the muscles twice. This is intended to relax them and boost circulation. The device is designed in particular for use on the upper back area (in particular the neck; **Picture B1 + B2**). But it can also be used to massage abdominal, buttock and leg muscles (**Picture B3–B6**).

Operating the Device

1. Connect the mains adapter to the connecting socket (**8**) of the device.
2. Insert the mains plug into a plug socket that is easily accessible.
3. Place the massage surface (**1**) of the device onto the part of the body which is to be massaged. For a shoulder or back massage, place the straps of the device around your shoulders and fix them in

place by putting your arms through the loops (7) (**Picture B7**) and thus holding the device in the correct position.

4. Start the massage by pressing the  button (3). Use the control panel (2) to select the function you want (see the "At A Glance" chapter).
5. After 15 minutes the device automatically switches off. To switch it off sooner, press the  button.
6. Pull the mains adapter out of the plug socket and disconnect the mains cable from the device.
7. Allow the device to cool down if the heat function was activated before it is stored away or cleaned (see the "Cleaning and Storage" chapter).

CLEANING AND STORAGE

Please Note!

- The device may not be connected to the mains for cleaning!
- Do not use any caustic or abrasive cleaning agents for cleaning. These can damage the surface.

- Wipe down the device just with a damp cloth if necessary. Then wipe it down fully with a dry cloth.
- Do not wind the mains cable around the device or the mains adapter.
- Keep the device in its respective packaging in a dry place which is protected from sunlight and is not accessible to children and animals.

TROUBLESHOOTING

If the device does not work properly, first check whether you are able to rectify the problem yourself. If the problem cannot be solved with the steps below, contact the customer service department.

Do not attempt to repair an electrical device yourself!

The device does not work.

The  button has not been pressed properly.

- ▶ Press the .

The mains adapter is not inserted correctly in the plug socket.

- ▶ Correct the fit of the mains adapter.

The device is not connected correctly to the mains adapter.

- ▶ Correct the connection.

There is too much weight on the device.

- ▶ Reduce weight on the device.

No heat is produced.

The red light and heat are not switched on.

- ▶ Press the  button.

The device gets too warm.

The device was in operation several times in succession for 15 minutes / too long in each case.

- ▶ Pull the mains adapter out of the plug socket and allow the device to cool down.

The device makes noises.

Operating noises are a normal consequence of the mechanics and not a defect.

Smoke emerges from the device.

The device is defective.

- ▶ Pull the mains adapter out of the plug socket immediately and contact the customer service department

TECHNICAL DATA

Article number: 14323
 ID operating instructions: Z 14323 M DS V1 0124 as

Device

Model number: C37-MP085
 Voltage supply: 12 V 
 Power: 24 W
 Timer: 15 minutes

Mains adapter

Manufactured for DS Produkte GmbH
 HRB 8937 HL
 Stormarnring 14, 22145 Stapelfeld, Germany
 Model number: LX-065-E120200
 Input: 100–240 V ~ 50/60 Hz; max. 0.6 A
 Output: 12 V DC  2 A; 24 W
 Protection class: II
 Ø efficiency in operation: 87.09 %
 Efficiency on low load (10 %): 81.82 %
 Power consumption on no-load: ≤ 0.10 W

DISPOSAL

Dispose of the packaging material in an environmentally friendly manner so that it can be recycled.

The adjacent symbol (crossed-out bin underlined) means that used devices should not be disposed of in the domestic waste, but should instead be taken to special collection and return systems.

Owners of used devices from private households can take them free of charge to the collection points of the public waste disposal authorities or the disposal points which are set up by manufacturers and distributors under the Electrical Equipment Act so that they can be disposed of in an environmentally friendly way and valuable raw materials can be recovered. If they are not disposed of properly, toxic ingredients may escape into the environment and cause harm to the health of humans, animals and plants. Businesses which market electrical and electronic devices are also obliged to take back used devices.

Waste prevention makes an even more valuable contribution to protecting the environment. This means if possible, as well as continuing to use or repair a device yourself, passing it on to another user also represents an ecologically beneficial alternative to disposal.

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All rights reserved.
Gefabriceerd voor DS Produkte GmbH